


# Compton Effect

Author: Bartu Yaman  
Lab Partner: Riza Serkan Kaskan  
Experiment Date: 25.03.2022  
Submission Date: 01.04.2022

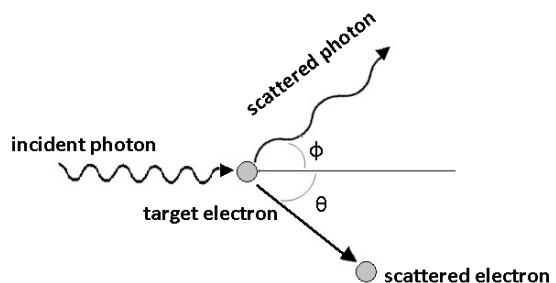

---

## Contents

<b>1</b>	<b>Introduction</b>	<b>4</b>
1.1	Goal . . . . .	4
1.2	Theory . . . . .	4
1.2.0.1	Thomson Scattering - Classical Electromagnetic Theory . . . . .	4
1.2.0.2	Compton Scattering - Particle Properties Of Light . . . . .	5
1.2.0.3	Bridging The Gap . . . . .	7
<b>2</b>	<b>Experimental Details</b>	<b>8</b>
2.1	Equipment . . . . .	9
2.1.1	Multi Channel Analyzer . . . . .	9
2.1.2	Radioactive Source . . . . .	9
2.1.2.1	CS - 137, 37 kBq . . . . .	9
2.1.2.2	CS - 137, 18.5 MBq . . . . .	10
2.1.3	Gamma Detector . . . . .	10
2.1.3.1	Scintillator . . . . .	10
2.1.3.2	Photomultiplier Tube . . . . .	10
2.1.3.3	Operating Unit . . . . .	11
2.1.4	Iron Rod . . . . .	11
2.1.5	Lead Block . . . . .	11
2.2	Procedure . . . . .	12
2.2.1	Calibration Procedure . . . . .	12
2.2.2	Experiment Procedure . . . . .	12
<b>3</b>	<b>Measurement And Data Analysis</b>	<b>14</b>
3.1	Presentation Of Data . . . . .	14
3.1.0.1	Calibration . . . . .	14
3.1.0.2	Main . . . . .	14
3.2	Calculations . . . . .	16
<b>4</b>	<b>Results And Discussion</b>	<b>17</b>
4.1	Discussion Of Results . . . . .	17
4.2	Discussion Of Errors . . . . .	17
4.2.1	Main Sources Of Errors . . . . .	17
4.2.1.1	Angle Arrangement . . . . .	17

4.2.1.2	Shielding	17
4.2.2	Error Mitigation	17

---

# 1 Introduction

## 1.1 Goal

The goal of this experiment is to observe the particle properties of light and to analyze the experimental observations to validate the wave-particle duality, bridging the gap between classical and quantum physics.

## 1.2 Theory

### 1.2.0.1 Thomson Scattering - Classical Electromagnetic Theory

Thomson scattering is an example of elastic scattering of electromagnetic radiation. The electric and magnetic fields components of an electromagnetic field exerts a Lorentz force on a particle, which sets it into motion that is periodic in time. Since the particle is now accelerated, it emits radiation.

The electric field component of a linearly polarized, monochromatic, plane wave incident on a particle of charge  $\mathbf{q}$  can be written as

$$E = eE_0 \exp(i(k.r - wt)) \quad (1)$$

The particle undergoes oscillations with small amplitude, therefore its velocity is assumed to be non-relativistic, enabling us to neglect the magnetic component of the Lorentz force. Then the EOM of the charged particle can be written as

$$F = qE = m\ddot{s} \quad (2)$$

where  $\mathbf{s}$  is the displacement from the origin.

Then the time averaged power radiated per unit solid angle becomes

$$\frac{dP}{d\Omega} = \frac{q^2 \langle \ddot{s}^2 \rangle}{16\pi^2\epsilon_0 c^3} \sin^2\theta \quad (3)$$


and

$$\langle \ddot{s}^2 \rangle = \frac{q^2}{m^2} \langle E^2 \rangle = \frac{q^2 E_0^2}{2m^2} \quad (4)$$

Which essentially means that the oscillating particle acts as a short antenna.

The resulting radiation will be perpendicular to the acceleration of the particle and be polarized along its direction of motion.

Figure 1: Thomson Scattering Diagram

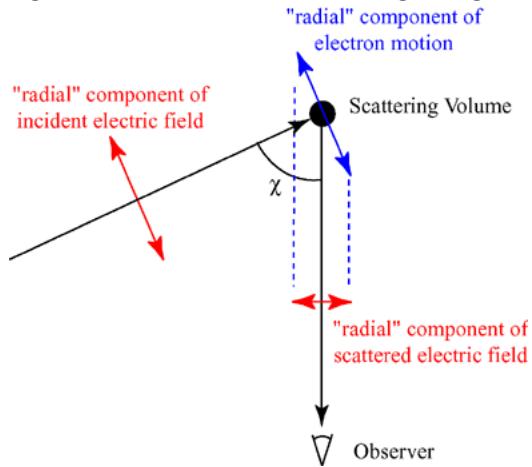

### 1.2.0.2 Compton Scattering - Particle Properties Of Light

Compton scattering is an example of inelastic scattering of electromagnetic radiation. Collision between a photon  $\gamma$  with wavelength  $\lambda$  and an electron  $e$  that is at rest causes the electron to recoil, which in turn results in a new photon  $\gamma'$  with wavelength  $\lambda'$  at angle  $\theta$  from the incident photon's path. Using the conservation of energy, we have

$$E_\gamma + E_e = E_{\gamma'} + E_{e'}. \quad (5)$$

Compton postulated that photons carry momentum and used the conservation of momentum as well

$$p_\gamma + p_e = p_{\gamma'} + p_{e'} \quad (6)$$


Knowing that,

$$E_\gamma = hf, E_{\gamma'} = hf' \quad (7)$$

and that,

$$E_e = m_e c^2, E_{e'} = \sqrt{(p_{e'} c)^2 + (m_e c^2)^2} \quad (8)$$

Figure 2: Compton Scattering Diagram

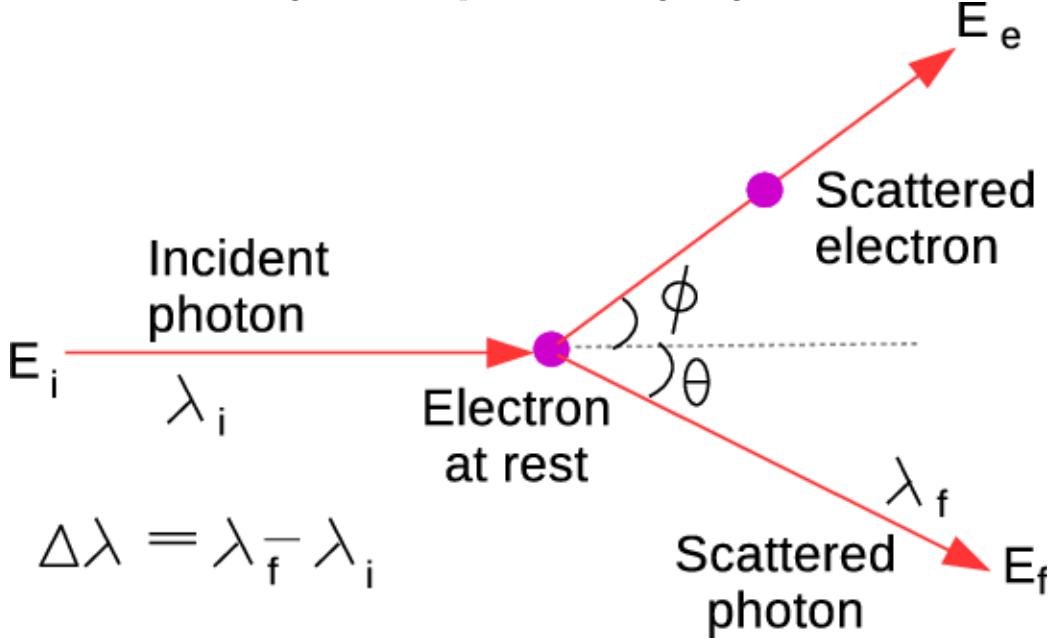

We expect the electron to be accelerated to a relativistic speed after the collision, thus we used the relativistic energy and momentum relation. Substituting these variables in the conservation relations, we get

$$hf + m_e c^2 = hf' + \sqrt{(p_{e'} c)^2 + (m_e c^2)^2} \text{ and } p_{e'}^2 c^2 = (hf - hf' + m_e c^2)^2 - m_e^2 c^4. \quad (9)$$

Solving the equation for the conservation of momentum,

$$\begin{aligned} p_{e'}^2 c^2 &= p_\gamma^2 c^2 + p_{\gamma'}^2 c^2 - 2c^2 p_\gamma p_{\gamma'} \cos\theta, \\ p_{e'}^2 c^2 &= (hf)^2 + (hf')^2 - 2(hf)(hf') \cos\theta, \\ 2hf m_e c^2 - 2hf' m_e c^2 &= 2h^2 f f' (1 - \cos\theta), \\ \frac{c}{f'} - \frac{c}{f} &= \frac{h}{mc} (1 - \cos\theta). \end{aligned}$$

And finally, after rigorous derivation, we arrive at the famous Compton Shift expression

$$\Delta\lambda = \lambda' - \lambda = \frac{h}{m_e c} (1 - \cos\theta). \quad (10)$$

Where  $\lambda$  is the wavelength of the incident photon,  $\lambda'$  is the wavelength of the scattered photon,  $m_e$  is the rest mass of the electron,  $\theta$  is the scattering angle and the expression  $\frac{h}{m_e c}$  is the Compton wavelength.

### 1.2.0.3 Bridging The Gap

As it is clear from the derivations, Thomson scattering is the low-energy limit of Compton scattering, where the photon energy is much smaller than the resting energy of the electron  $\mu \ll mc^2/h$ . Since classical electrodynamics fails to explain the wavelength shift between the incident and scattered photons, realization of this phenomena had important implications on physicists' understanding of light since light clearly behaves as if it consists of particles in this experiment.


## 2 Experimental Details

The experiment was carried out with utmost care, following the radiation safety guidelines. The surroundings of the path of the photons radiated from the gamma radiation sources were carefully obstructed and the experiment area was checked with a Geiger counter before anyone interacted with the setup.

The procedure consisted of two parts that will be explained individually later in this section. The experiment took approximately two hours including the two procedures, discussion about the experiment with the TA and necessary calculations.

The experiment setup used can be seen in the figure below, the equipment will be explained individually.

Figure 3: Experiment Setup

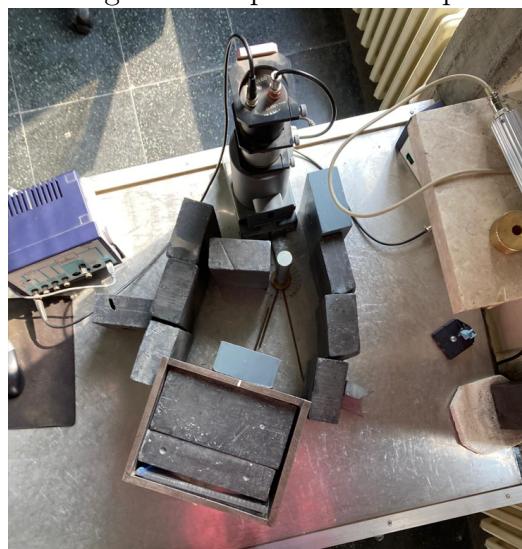


## 2.1 Equipment

### 2.1.1 Multi Channel Analyzer

A multichannel analyzer is an instrument that uses a fast analog to digital converter to analyze incoming input pulses and is a popular instrument used in gamma spectroscopy. The data taken from the MCA will be presented in the next section.

Figure 4: The Multi Channel Analyzer Used In The Experiment

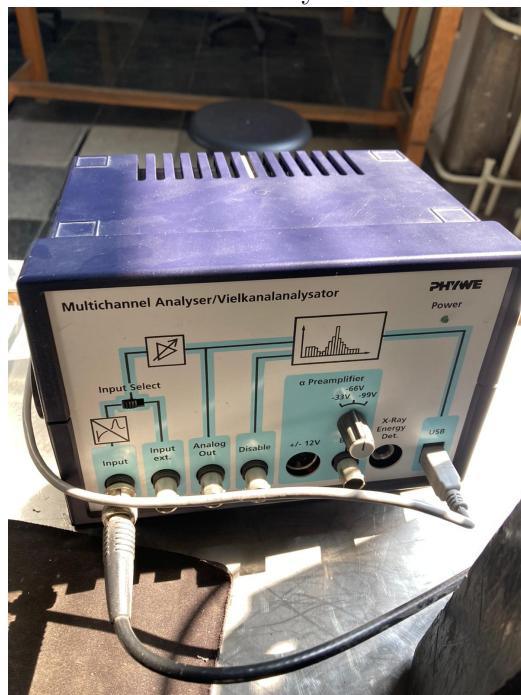

### 2.1.2 Radioactive Source

Two radioactive sources were used during the experiment.

#### 2.1.2.1 CS - 137, 37 kBq

This source was encased in a small lead block and used during the calibration process. Its small size allowed it to be stored safely without much effort. When the calibration process was over, this source was carefully removed from the setup and left adjoining the wall.

### 2.1.2.2 CS - 137, 18.5 MBq

This source was encased in a large lead block, had wheels attached to it and used during the experiment procedure. A thick lead block was present in front of this source during the calibration process for safety reasons, the block was only removed briefly for the experiment procedure. The source is encased in the large lead block which can be seen in Figure 3.


### 2.1.3 Gamma Detector

The gamma detector used in this experiment is a scintillation counter, which consists of a scintillator crystal and a photomultiplier tube connected to it.

#### 2.1.3.1 Scintillator

Scintillators absorb incoming particles (or ionizing radiation, ie. gamma radiation used in this experiment) to emit the absorbed energy as light. Scintillators themselves work on the principle of Compton effect, which is a nice correlation.

Figure 5: A Scintillator Crystal

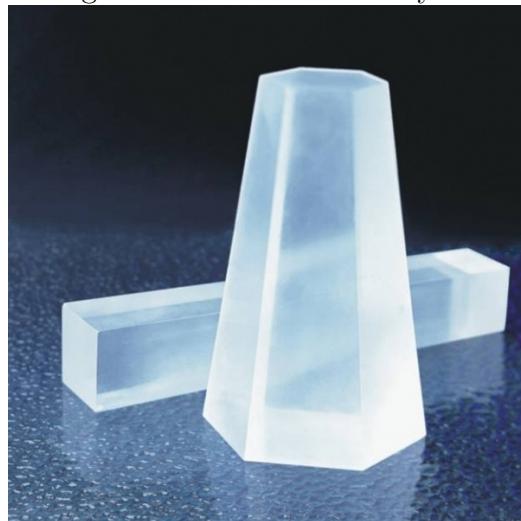


#### 2.1.3.2 Photomultiplier Tube

A photomultiplier tube allows the incoming photons to hit the photocathode material inside, allowing electrons to eject its surface because of photoelectric

effect. The ejected electrons are then directed over to the electron multiplier where secondary emission allows the electrons to be multiplied.

Figure 6: Photomultiplier Tube Schematic

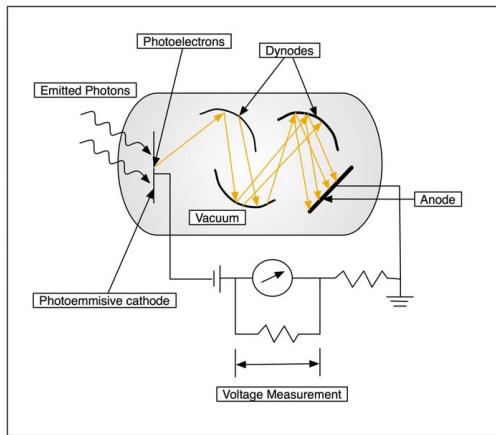

### 2.1.3.3 Operating Unit

The operating unit of the detector is responsible for supplying the appropriate voltage in a stable way to the photomultiplier tube. It operates between 600 Volts and 1100 Volts.

During the experiment, the voltage knob was at the scale 8.00, which corresponds to 1000 Volts.

### 2.1.4 Iron Rod

An iron rod is used as an electron source in this experiment. The photons radiated by the gamma source hit the iron rod and scatter photons to be detected by the detector.

### 2.1.5 Lead Block

Several lead blocks are used in the experiment for safety reasons. Everywhere except the direct path of the incident radiation is obstructed by lead blocks, which in turn has consequences in the measured result which will be discussed in the discussion section.

## 2.2 Procedure

### 2.2.1 Calibration Procedure

The calibration process was performed as follows:

- The MCA was turned on and its proprietary software was set up to make measurements.
- The voltage knob of the operating unit of the detector was adjusted to the scale value of 8.00 to supply the necessary voltage to the detector.
- Cs - 137, 37kBq, which is the source encased in the smaller lead block, was carefully inserted in front of the detector, adjoining it.
- The calibration process was started in the software.
- Waited around 15 - 20 minutes for the measurements.
- Analyzed the resulting spectra with the TA.
- Identified the two sharpest peaks corresponding to two different types of radiation from the source and marked them.
- Assigned known peak values for Cs - 137 to the marked peaks, completing the calibration process.
- The calibration was saved on the software and the radiation source was carefully removed.

### 2.2.2 Experiment Procedure

The experiment procedure was performed as follows:

- The software was set up to make measurements again.
- The iron rod, which acts as an electron source, was inserted in front of the larger gamma source and its angle was arranged with an old and worn down instrument.
- The direct path of the radiation was blocked by a lead block to ensure that we mostly measure the radiation due to Compton scattering.

- The lead cover of the Cs - 137, 18.5 MBq source was carefully removed to let radiation hit the iron rod.
- The measurements was started with the software.
- Waited around 15 - 20 minutes one again for the measurements.
- After we had enough data, the lead cover of the radiation source was carefully inserted in front of it once again and the measurement process was done.
- The resulting spectra was analyzed with the TA and had verbal discussion on it.
- Performed smoothing on the measured data and saved the plots for the laboratory report.

## 3 Measurement And Data Analysis


### 3.1 Presentation Of Data

Since we are only interested in the peak values for our calculations and the general shape of the spectra to draw conclusions, the data is taken as photographs from the laboratory computer.

#### 3.1.0.1 Calibration

The calibration spectra data is presented below, the significance of the peaks and their values will be discussed in the discussion section.


Figure 7: Data Taken During The Calibration Process

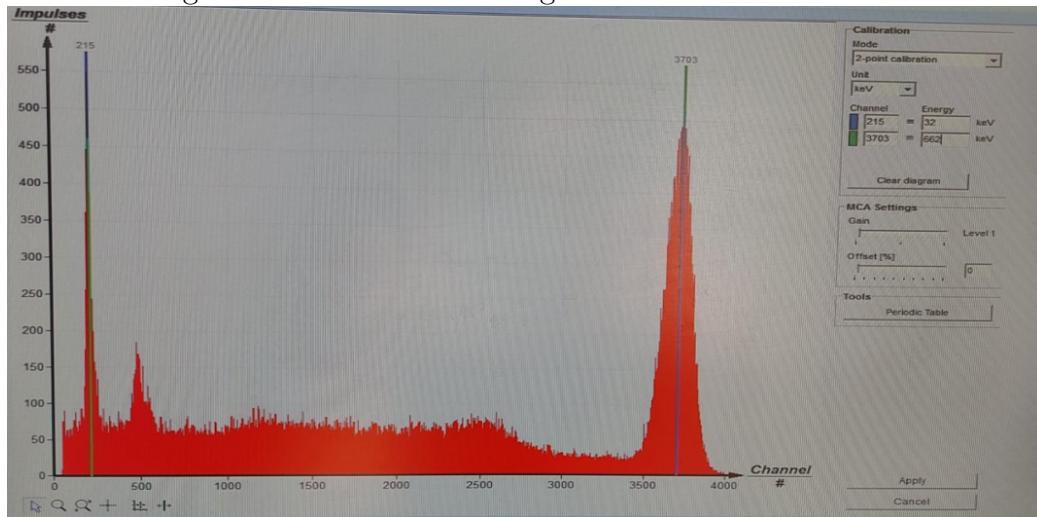

#### 3.1.0.2 Main


The raw spectra data from the experiment is presented below, this spectra was then smoothed for easier analysis. From the figure, we can see that the measurement was taken for 1331 seconds and there were 336247 total number of impulses.

Figure 8: Raw Data Taken During The Main Part

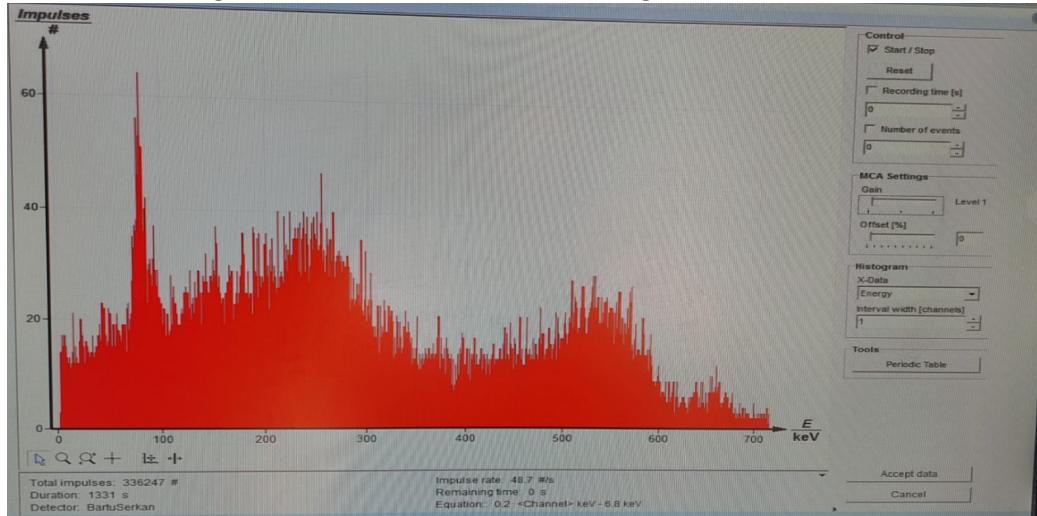

The smoothed spectra data from the experiment is presented below, the meanings of the three significant peaks will be discussed in the discussion section.

Figure 9: Smoothed Data Taken During The Main Part

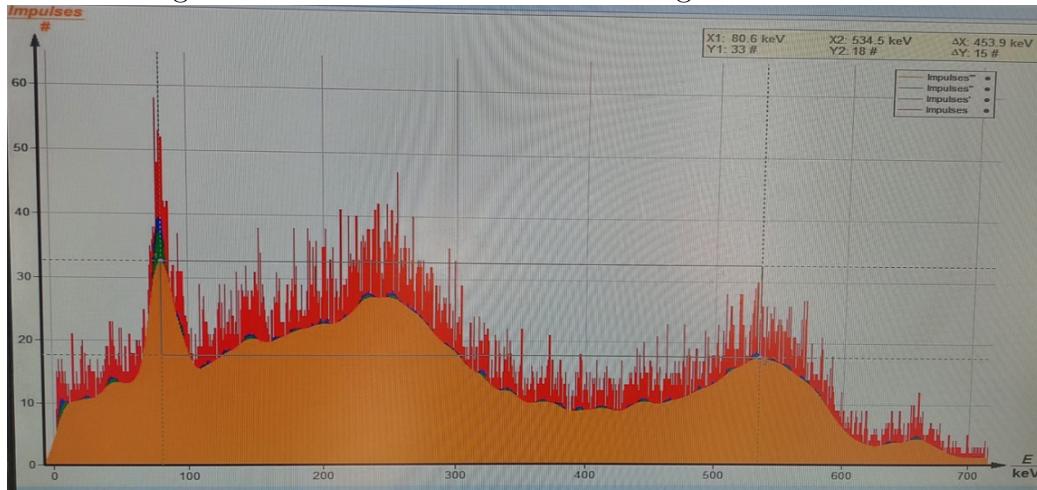

### 3.2 Calculations

For the proof of concept, we needed to calculate the Compton wavelength which we initially knew the value of. The calculations were carried out as follows:

$$\text{incident photon wavelength, } \lambda = 1.873 * 10^{-12} m$$

$$\text{incident photon energy, } E = 661.6 keV$$

$$\text{scattering angle, } \phi = 30^\circ$$

$$\text{rest mass of electron, } m_0 = 0.510998950 MeV/c^2$$

$$\frac{1}{E'} - \frac{1}{E} = \frac{1}{m_0 c^2} (1 - \cos\phi)$$

$$E' = 563.82 keV$$

$$E' = hc/\lambda'$$

$$\lambda' = 0.0219928 \text{\AA}$$

$$\lambda = 0.0187311 \text{\AA}$$

$$\lambda_c = \frac{\lambda' - \lambda}{(1 - \cos\theta)} = 0.0243457 \text{\AA}, \text{ the Compton wavelength}$$

## 4 Results And Discussion

### 4.1 Discussion Of Results

### 4.2 Discussion Of Errors

Although some error was present in the measurement, the results were clear enough to draw concrete conclusions. Especially after the smoothing done by the software of the MCA, errors were less prominent.

#### 4.2.1 Main Sources Of Errors

The main sources of error in this experiment are:

##### 4.2.1.1 Angle Arrangement

##### 4.2.1.2 Shielding

#### 4.2.2 Error Mitigation

## References

- [1] IBM. Ibm quantum composer docs, 2021.
- [2] Michael A. Nielsen and Isaac L. Chuang. *Quantum Computation and Quantum Information*. Cambridge University Press, 2010.